

Gypsums



Prestige Dental Products, Inc

Prestigedentalproducts.com

Tel: 877-772-3888 Fax: 714-632-8821

Gypsums

Disposable Model Base Formers (Dentsply)




Saves clean up time. Assures creation of consistent models. Durable plastic will not deform during model creation. Flexible plastic for easy separation.

60/Pkg,
20 each small, medium and large 34-11111

120/Pkg,
40 each small, medium and large 34-11112

Disposable Model Base Holders (PlasDent)

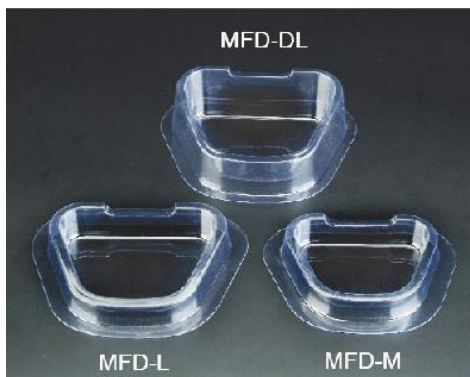

Disposable Large/Deep, Clear, 60
pcs/box 34-11132

Disposable Large, Clear,
60pcs/box 34-11131

Disposable Medium Clear, 60pcs/
box 34-11130

**Disposable Asst Sizes 20 of each
size,** 60pcs/box 34-11129

**Disposable Asst Sizes 40 of each
size,** 120pcs/box 34-11128

Gleco Plaster Trap (IndiGreen)

Trap System Invented by a plumber, the Gleco Trap brings revolution to the age-old chore of sink trap retrieval and plaster trap maintenance. No more leaky bags or bulky buckets!

Install the pre-assembled Gleco Trap less than one minute start to finish without tools. Simply turn side valve to drain off excess water, unscrew full bottle, and screw in a clean one.

- Side drain siphons out all water prior to removal, eliminating spillage
- Translucent bottle shows waste level.



Prestige Dental Products, Inc.

Prestigedentalproducts.com

Tel: 877-772-3888

- Adapts easily to existing plumbing — the Gleco Trap is the only plaster trap to meet US Uniform Plumbing Code!

System includes one trap assembly, two replacement bottles, all necessary 1 1/2" pipe components and simple instructions for a five-minute installation

64 oz. Gleco Trap System 34-87931

19 oz. Gleco Trap System 34-87932

**19 oz. replacement bottle for
Gleco, 6/Box** 34-87933

Gleco Casket, 2/Pkg 34-87934

Gleco Trap HV System (IndiGreen)



Item #10012 64 oz. Gleco Trap System
Item #10010 19 oz. Gleco Trap System

Item #11010 3.5-Gallon Gleco Trap HV Kit
Item #11011 5-Gallon Gleco Trap HV Kit

From the inventor of the original Gleco Trap, the Gleco Trap HV Sanitary High-Volume Kit brings the same clean, dry and simple maintenance to full-size plaster traps. Never again will you soak your cabinets and yourself with foul-smelling, contaminated wastewater when you change your plaster trap container.

The Gleco Trap HV features a specially-engineered air-tight lid with threaded inlet and outlet extensions that reserve an air pocket inside. When the trap is full, simply remove the air-release plug to drain all inlet and outlet plumbing into the bucket before disconnecting. Its patented design prevents the mess and contamination that occurs with the removal of traditional plaster traps.

Instead of leak-prone flexible hoses found on other traps, the Gleco Trap HV uses standard 1 1/2" threaded inlet and outlet ports to ensure a positively leak-proof assembly. Adapts easily to existing plumbing.

Each kit includes a heavy-duty container, all 1 1/2" PVC plumbing components to fit most sink configurations and instructions.

1 gallon, Refill (4) 34-87945

32 oz. Refill (12) 34-87946

64 oz. Refill (6) 34-87947

GT-19 w/19oz. Bottle 34-87948

GT-2 w/64oz. Bottle 34-87950

Replacement Gasket Set 34-87950



Prestige Dental Products, Inc.
 Prestigedentalproducts.com
 Tel: 877-772-3888

Green Rubber Base Formers for Plastic Articulators



Q-Left	34-11117
Q-Right	34-11116
Small	34-11115
Medium	34-11113
Large	34-11114
X-Large	34-11118

Green Silicon Rubber Base Formers



Base formers made of silicon rubber (Left, Right, Small, Medium, Large). Simply pour stone, level off, then place the teeth cast on the form Easy to use and reusable

[TOP: Medium full arch; Large full arch]
 [BOTTOM: Small full arch; Right & Left Quadrant]

Set of 5/Box	34-11119
---------------------	----------

Magnum Orthodontic Stone (MDC)

Magnum® Ortho stone

Class II

Snow white
 Specially designed for orthodontic models
 Super smooth and creamy



CHARACTERISTICS
 • Class III super stone
 • Regular setting
 • Extra fine powder
 • Optimal hardness

ADVANTAGE
 • Smooth surface
 • Creamy smooth gypsum
 • High detail
 • Low setting expansion
 • Dimensionally stable

PHYSICAL PROPERTIES
 Setting expansion 0.15%
 Pressure resistance 8000 psi
 Setting time 10 ± 2 min.
 Water/Powder ratio 32 cc / 100 gr

Presentation	Code
Box 22 lbs / 10 kg	22OP02



Magnum Orthodontic Stone, 22 lbs.	34-11120
--	----------

Gypsums

Magnum Plaster White Stone (MDC)

Magnum Plaster

Class II



CHARACTERISTICS

- Class II plaster
- Regular setting
- Extra fine powder
- Optimal hardness

ADVANTAGE

- Smooth surface
- Creamy gypsum
- High detail
- Low setting expansion
- Dimensionally stable

PHYSICAL PROPERTIES

Setting expansion	0.20%
Pressure resistance	5000 psi
Setting time	6 ± 2 min.
Water/Powder ratio	24 cc / 100 gr



Presentation	Code
Box 22 lbs / 10 kg	20PW02

Creamy White Stone, 22 lbs.	34-11121
------------------------------------	----------

Magnum Yellow & Blue Stones (MDC)

Magnum®

Stone class III

Manufactured under ISO and ADA specifications
 Also ISO 9001 - 2008 standard.



CHARACTERISTICS

- Class III super stone
- Regular setting
- Extra fine powder
- Optimal hardness

ADVANTAGE

- Smooth surface
- Creamy smooth gypsum
- High detail
- Low setting expansion
- Dimensionally stable

PHYSICAL PROPERTIES

Setting expansion	0.15%
Compression resistance	390 kg/cm2
Setting time	10 ± 2 min.
Water/Powder ratio	24 cc / 100 gr




Presentation	Yellow	Blue
Box 22 lbs / 10 kg	26Y02	22B02

Yellow, 22 lbs.	34-11122
------------------------	----------

Blue, 22 lbs	34-11123
---------------------	----------

Model Base Formers (Plasdent)

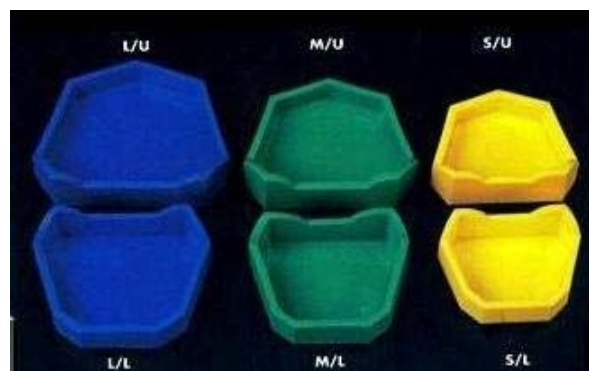


Large Blue, 3pcs upper & 3pcs lower/Box	34-11127
--	----------

Gypsums

Medium Green , 3pcs upper & 3pcs lower/Box	34-11126
Small Yellow , 3pcs upper & 3pcs lower/Box	34-11125
Asst Sizes , 6pcs/Box	34-11124

Model Carriers

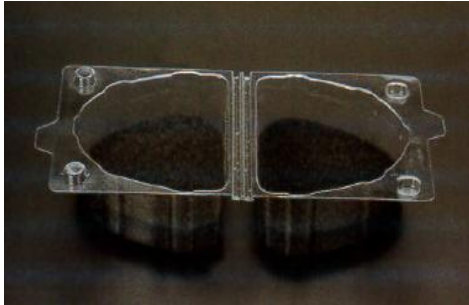

Carrier Set: 1 tray w/2 foams, 300 set/case	34-98041
Foams for model carriers, Charcoal, 1000/Pkg	34-98042

Trap-Eze Disposable Plaster Traps (Buffalo)

Economical trap systems designed to slow water discharge flow

- Requires minimum 9" clearance between bottom of sink tail pipe and top of trap bucket for easy installation without tools
- All trap buckets feature high-strength plastic construction and permanently-sealed lids
- 3 1/2 gallon size (12 1/2" diameter x 11" H) ideal for intermediate applications
- 5 gallon size (12 1/2" diameter x 14 1/2" H) best for larger offices and labs with higher volume collection needs (space permitting)
- Two flexible hoses (1 1/2" ID x 1 3/4" OD x 20" long) conveniently connect trap to sink tail pipe and sewer tail pipe
- Flexible adaptor fitting for sewer tail pipe included for improved drain seal
- Complete kit includes sealed trap bucket, 2 hoses, and sewer tail pipe adaptor fitting
- Refill includes trap and two plastic caps for sealing full trap
- Fast and easy installation and trap replacement
- Made in the USA



3 1/2 Gal. Trap-Eze Complete Kit (11" Height Sealed Bucket, Hoses)	34-11133
--	----------



Prestige Dental Products, Inc.
Prestigedentalproducts.com
Tel: 877-772-3888

3 1/2 Gal. Trap-Eze Refill (11" Height Sealed Bucket)	34-11134
5 Gal Trap-Eze Complete Kit (14.5" Height Sealed Bucket, Hoses)	34-11135
5 Gal. Trap-Eze Refill (14.5" Height Sealed Bucket)	34-11136
3 1/2 Gal. Trap-Eze SS Self-Sanitizing Trap Complete Kit New	34-11137
3 1/2 Gal. Trap-Eze SS Self Sanitizing Trap Complete Kit-New	34-11138
5 Gal Trap-Eze SS Self-Sanitizing Trap Complete Kit-New	34-11139
5 Gal Trap-Eze SS Self Sanitizing Trap Refill- New	34-11140
Hoses only for Trap-Eze, Set of 2 (18" Overall Length Each) Flexible	34-11141
Drain Fitting for Trap-Eze Exhaust Hose	34-11142
Trap-Eze BT Bottle Trap Complete Kit (1 1/2" Plumbing + One 64 oz. Bottle	34-11143
Trap-Eze BT Bottle Trap 64 oz. Re-fills w/Lids & Bibs, Pkg of 2	34-11144

www.PrestigeDentalProducts.com




**Check our
website for more
items and
variation**



Prestige Dental Products, Inc.
Prestigedentalproducts.com
Tel: 877-772-3888

Gypsums

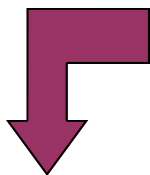

IndiGreen
INNOVATIONS

Brighter Ideas for Better Dentistry



Buffalo
Buffalo Dental Manufacturing Co., Inc.



Call Us To Place
Your Order Or Visit
Us Online at
www.prestigedentalproducts.com

Order Now



Call us for Quarterly Promotion

Prestigedentalproducts.com



Prestige Dental Products, Inc.
www.prestigedentalproducts.com

Order Form

Date:

Ordered By

Deliver To

☐ Same as Above

Business Name:

Address:

State/Province:

Zip/Postal Code:

Phone:

Fax:

Contact Name:

Business Name:

Address:

State/Province:

Zip/Postal Code:

Phone:

Fax:

Contact Name:

Item Number	Description	Quantity	Unit Price	Amount

Payment

☐ Credit Card

☐ Online Payment (Bank to Bank)

E-mail

☐ American Express

☐ MasterCard

☐ Visa

Sub Total

C.O.D. Charges

Freight Charges

Sales Tax

Total

☐ Check Payable To

Card Number:

Expiration Date:

Cardholder Name:

Dental License #

Year Licensed

Expiration Date

DEA #

Year Licensed

Expiration Date



Prestige Dental Products, Inc

3031 E. Coronado St., #D

Anaheim, Ca 92806

Tel:877-772-3888 Fax:714-632-8821

www.prestigedentalproducts.com

